



Shortcuts and on frame interactions



How do you activate shortcuts?

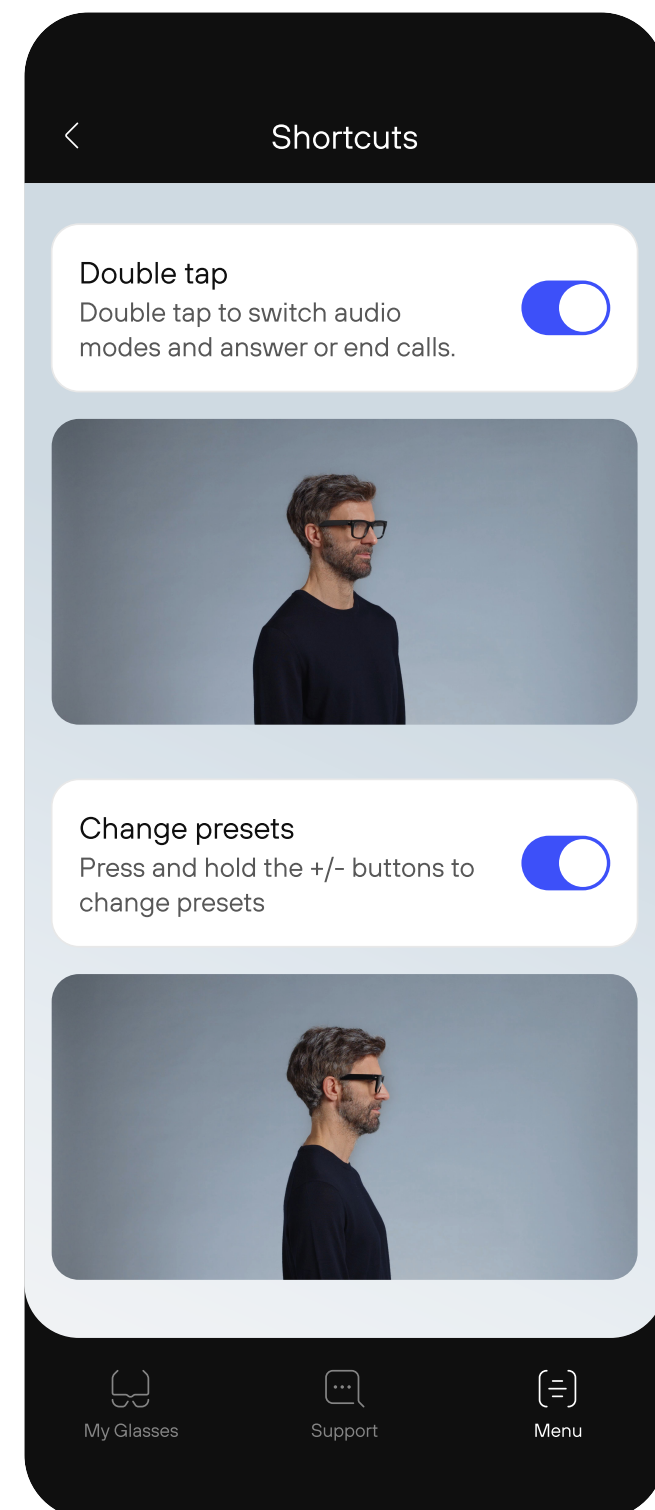

Go to **Menu > Shortcuts >** toggle the switch on "**Double tap**" and "**Change presets**".

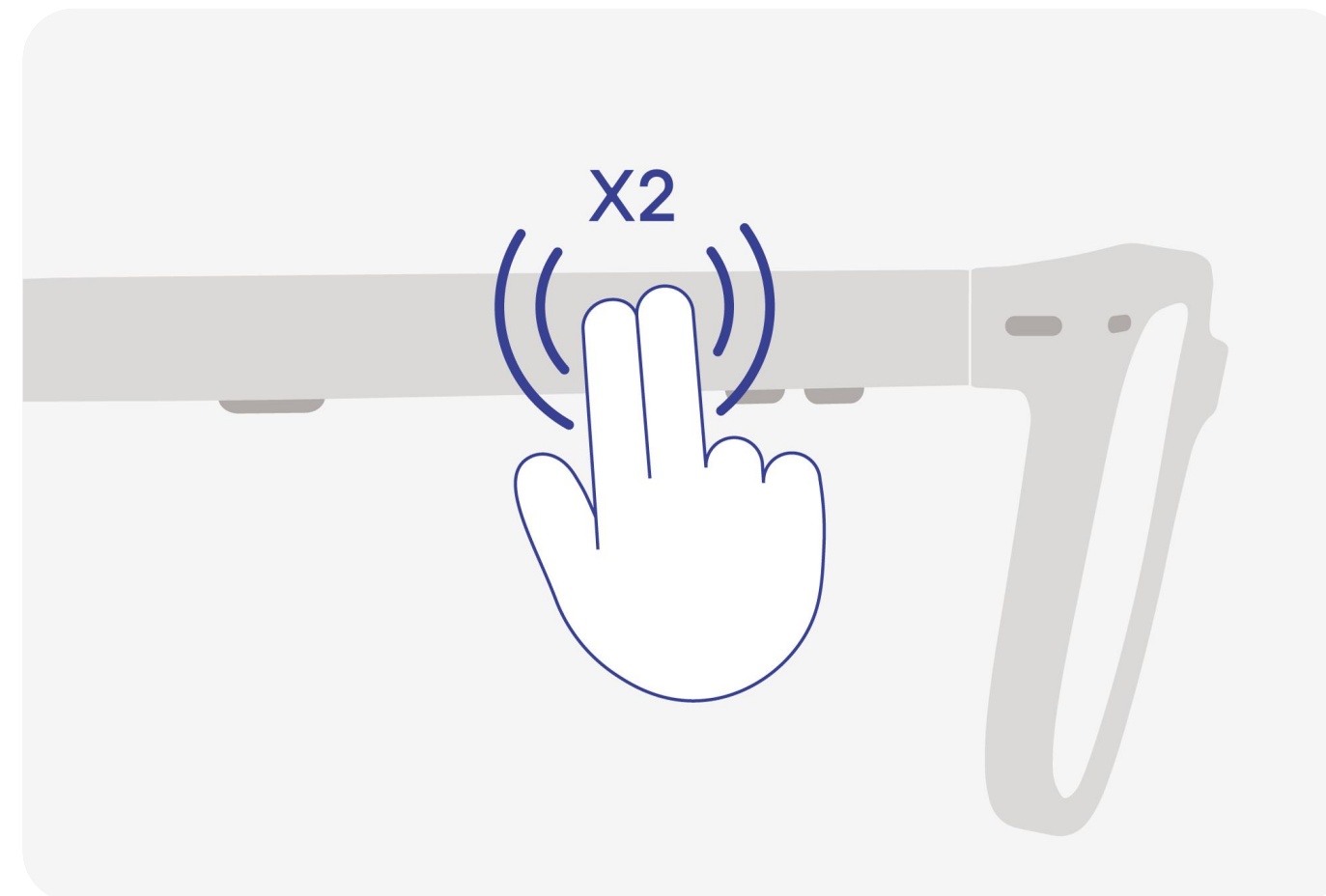




Shortcuts on the frame

Hearing Assistance (HA)


Switch Audio mode
Double tap (Opt-in).

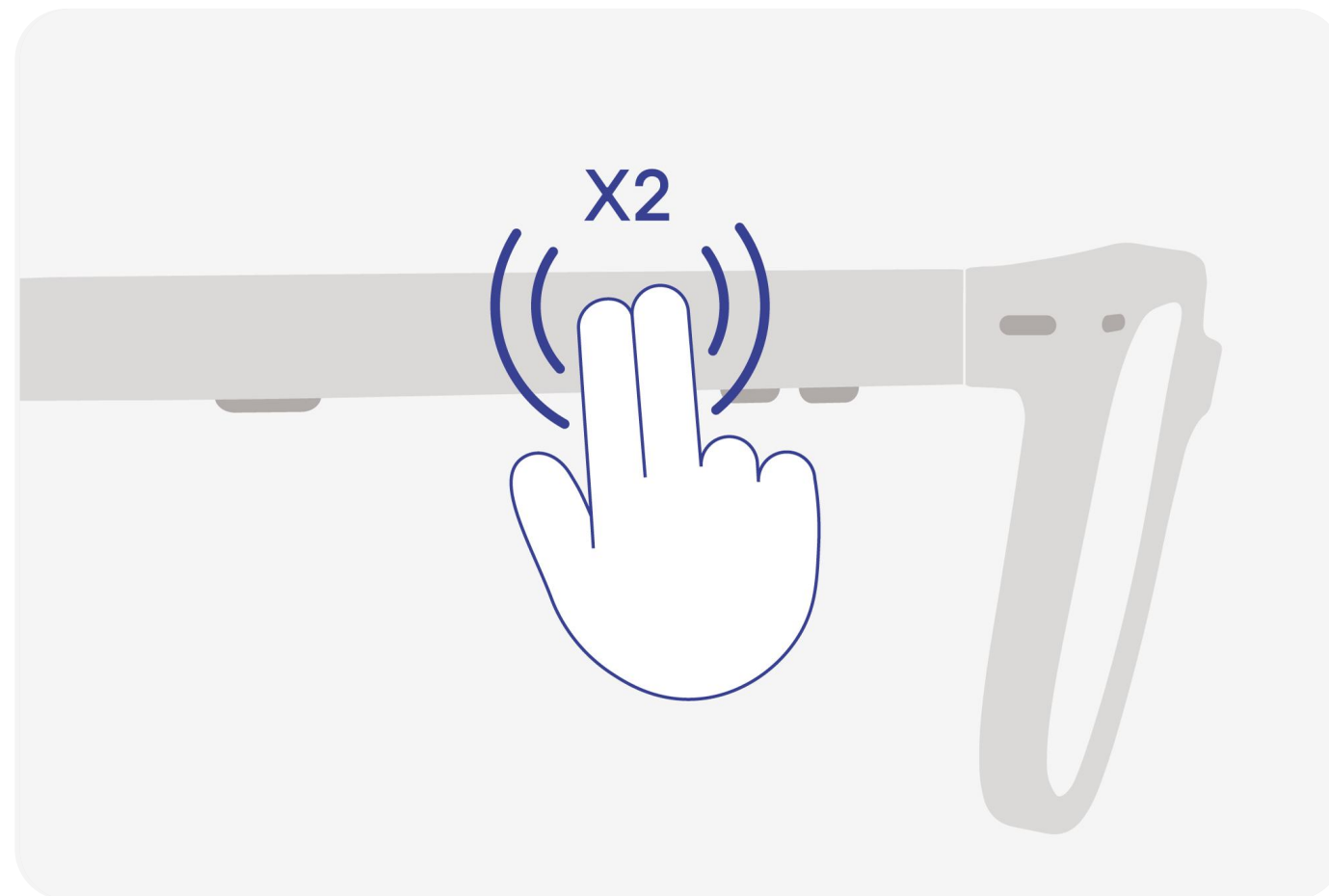

Preset up
Press (+) and hold 1.25 seconds (Opt-in).



Preset down
Press (-) and hold 1.25 seconds (Opt-in).

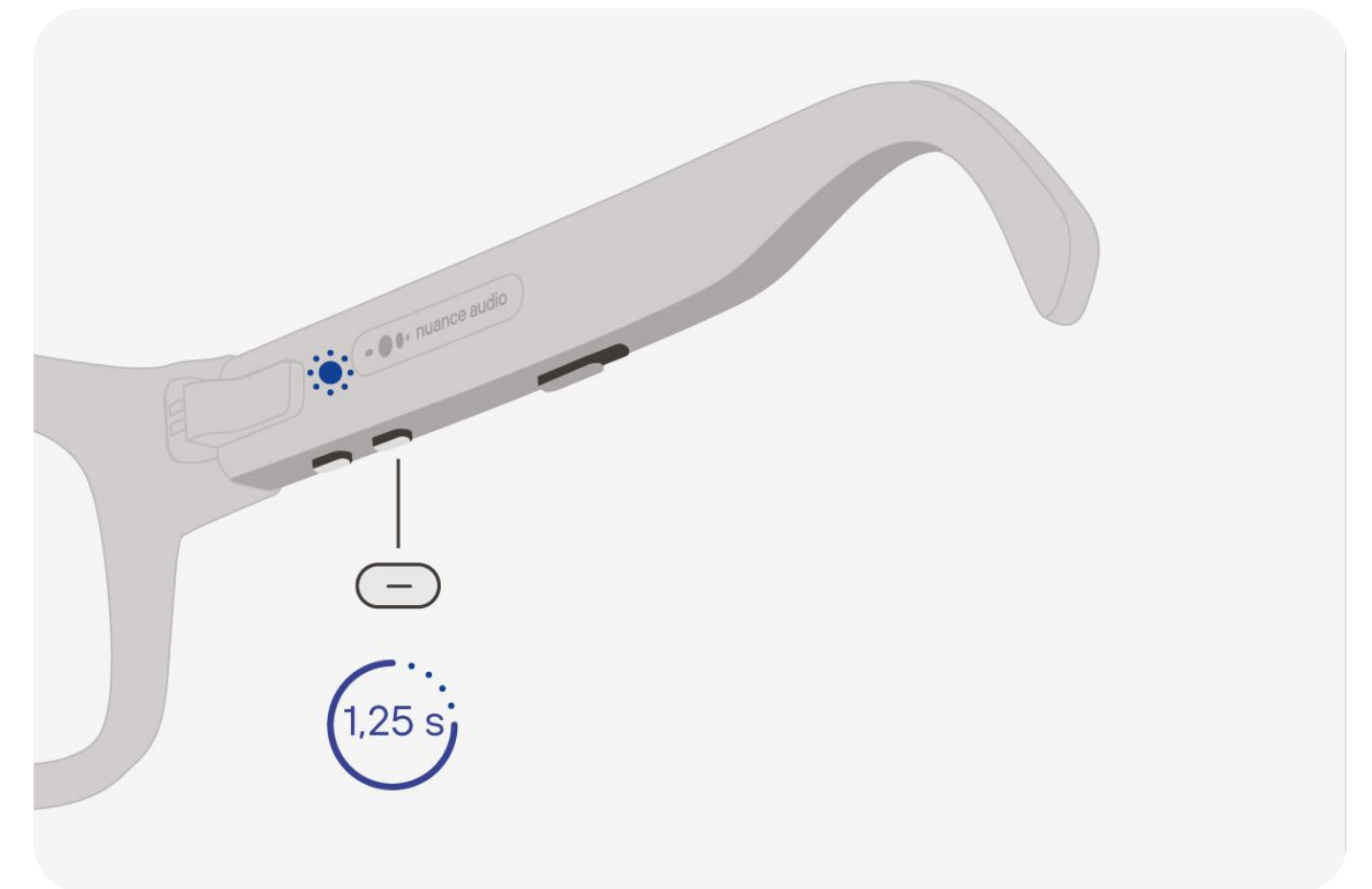





Shortcuts on the frame

Calls

Answer
Double tap (Opt-in).

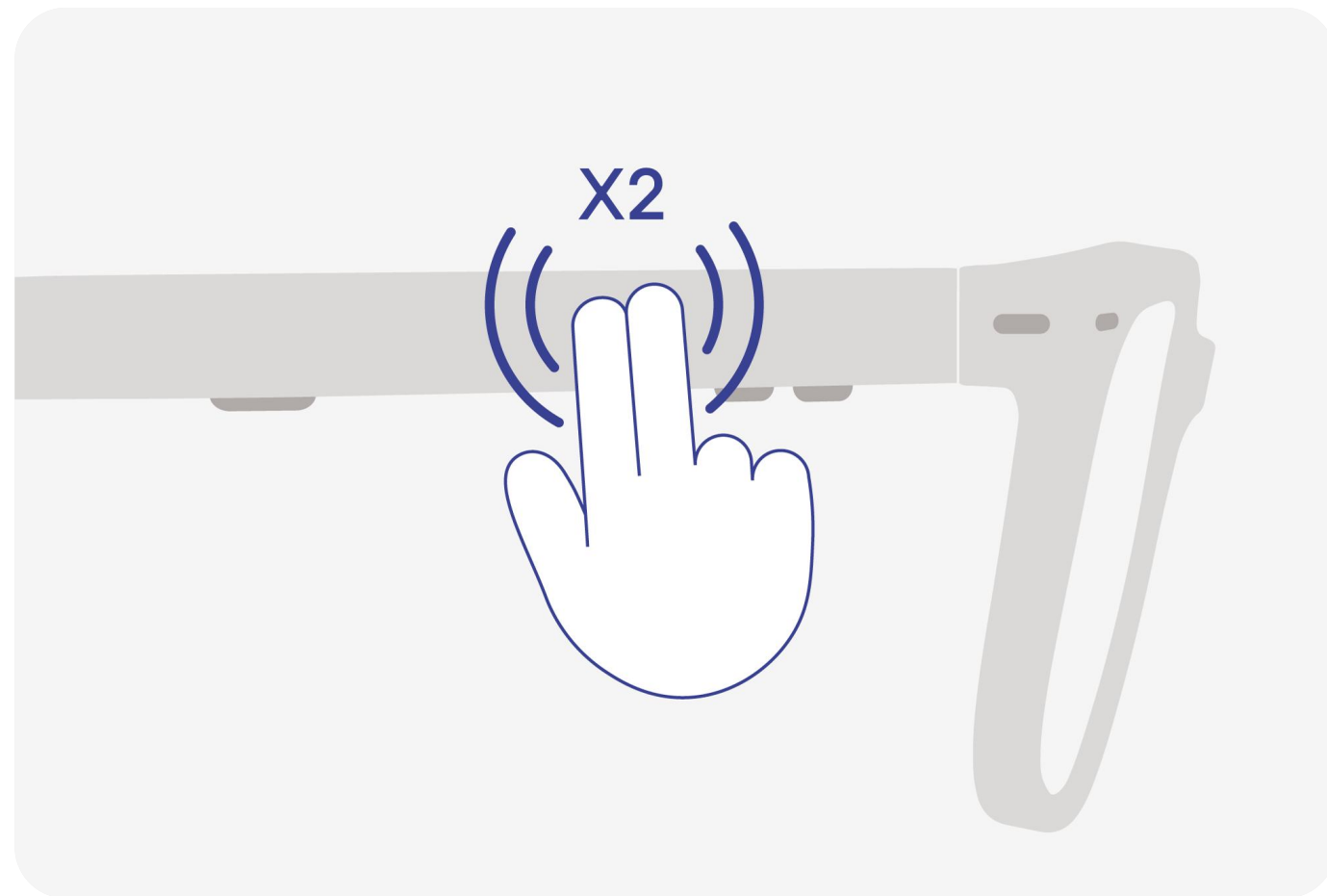

End
Double tap (Opt-in).

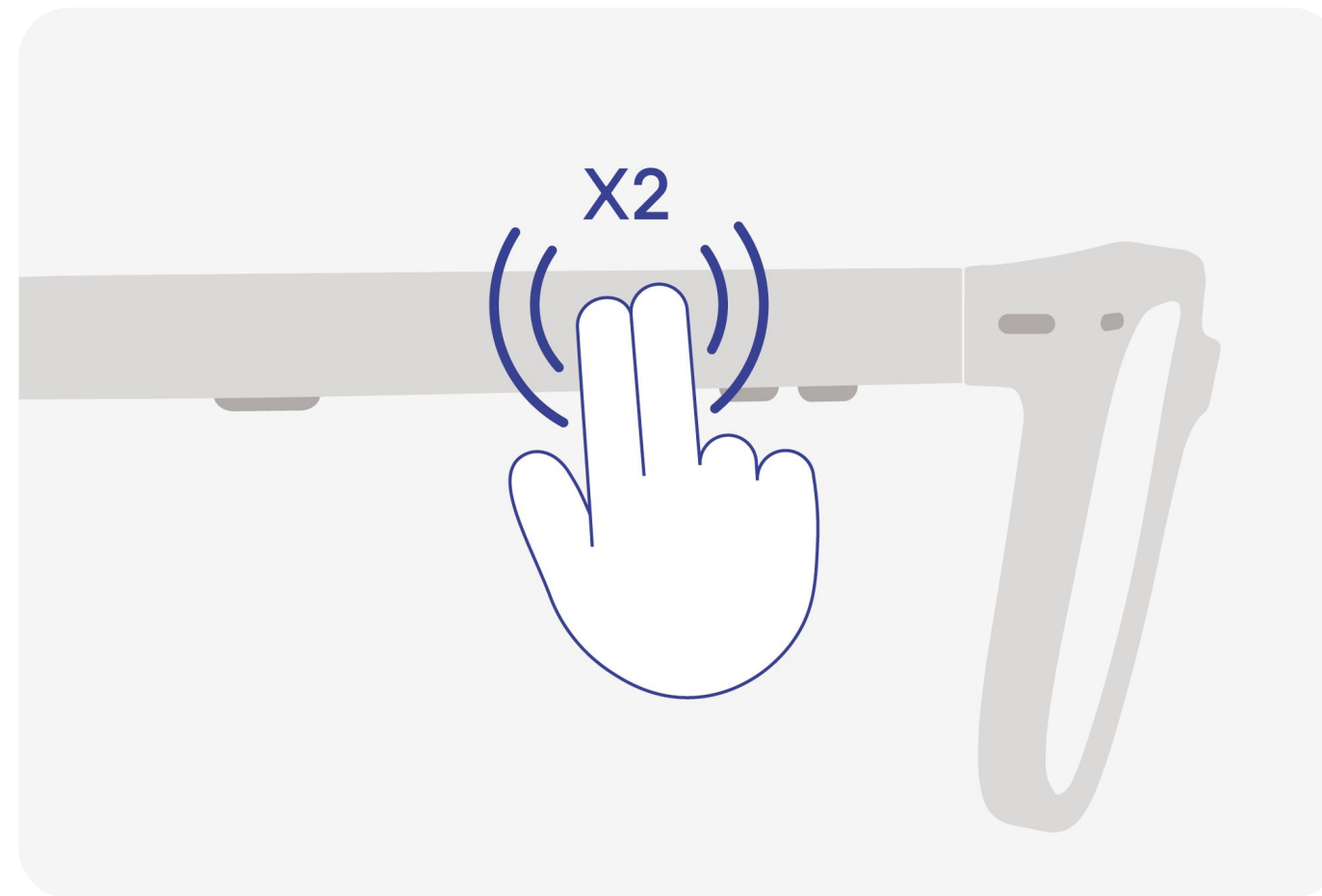





On frame interactions

General


Power on
Slide on

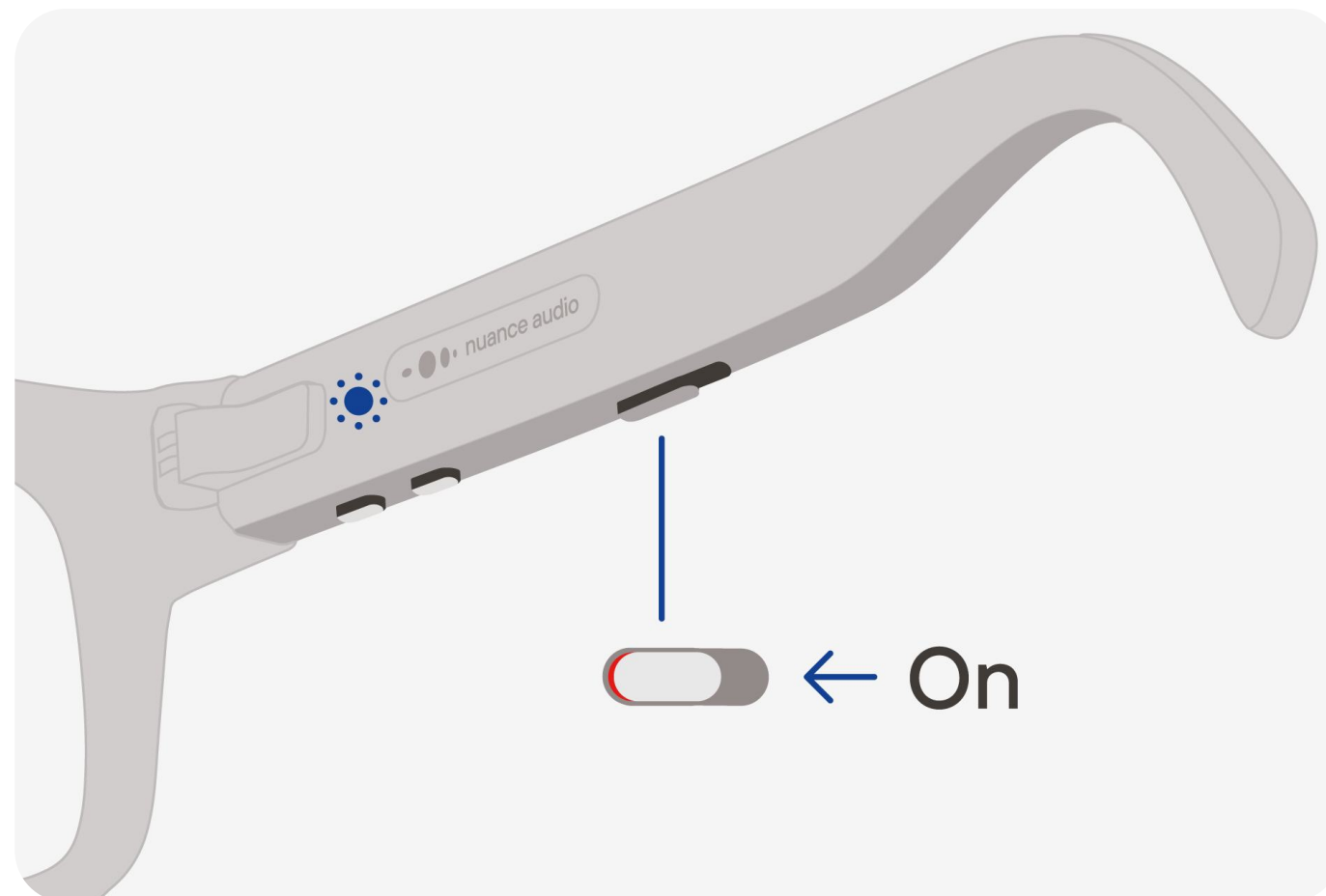

Power off
Slide off

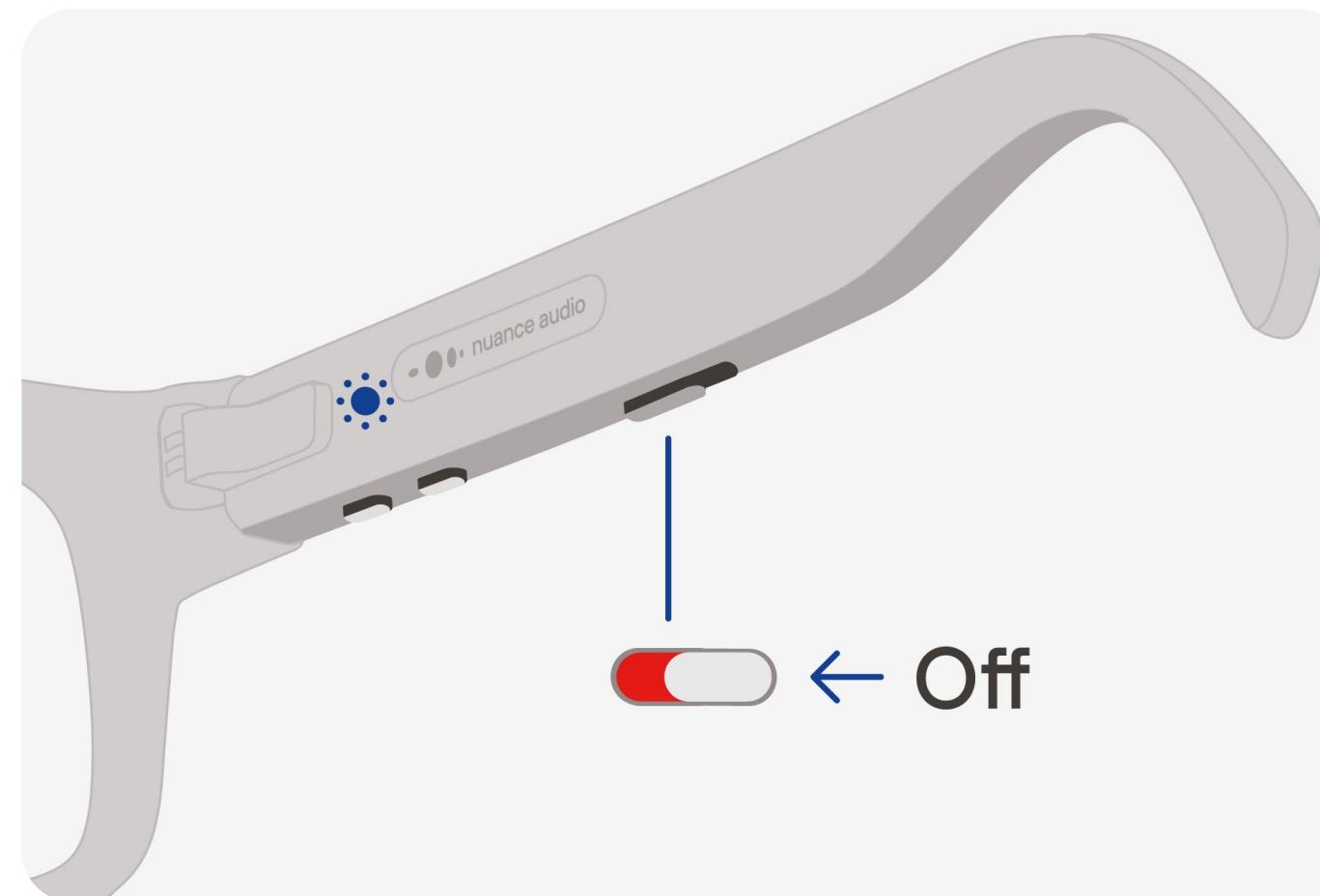

Force shutdown
Slide off then press (-) and hold 12 seconds.

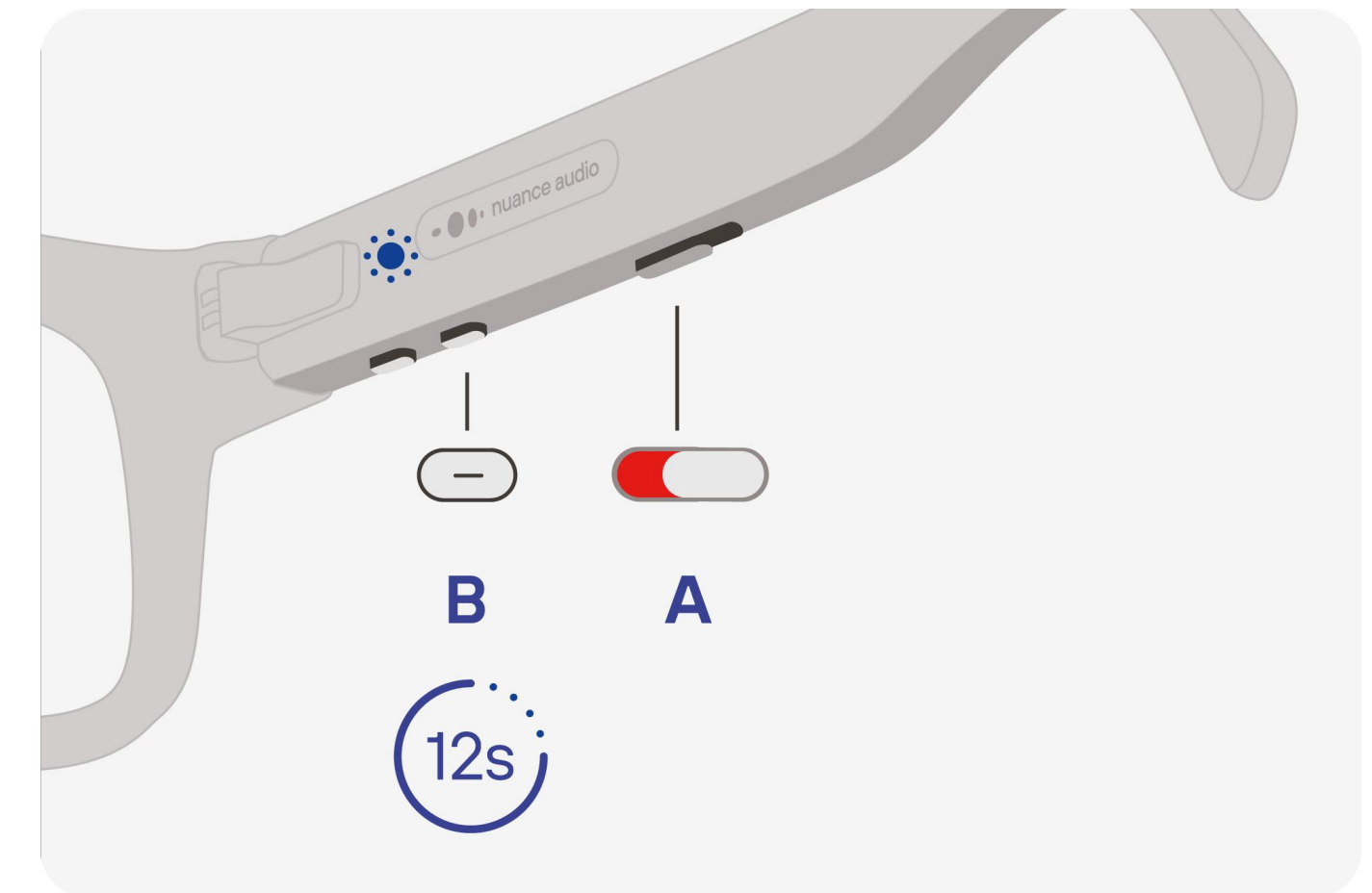




On frame interactions

General

1st Pair

Slide on

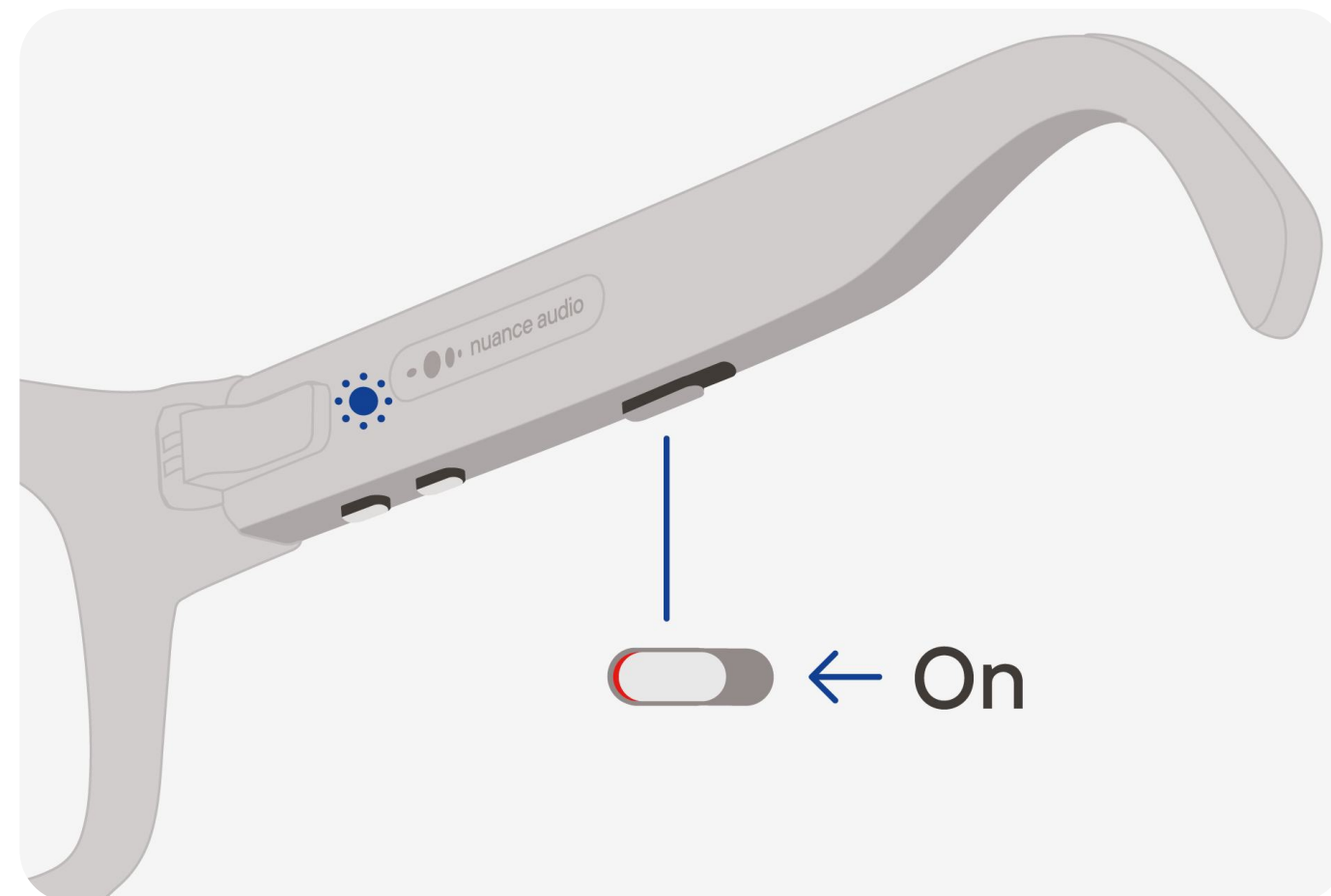

2nd Pair

Press (+) or (-) and hold for 8 seconds.






On frame interactions

General


Exit

Slide off

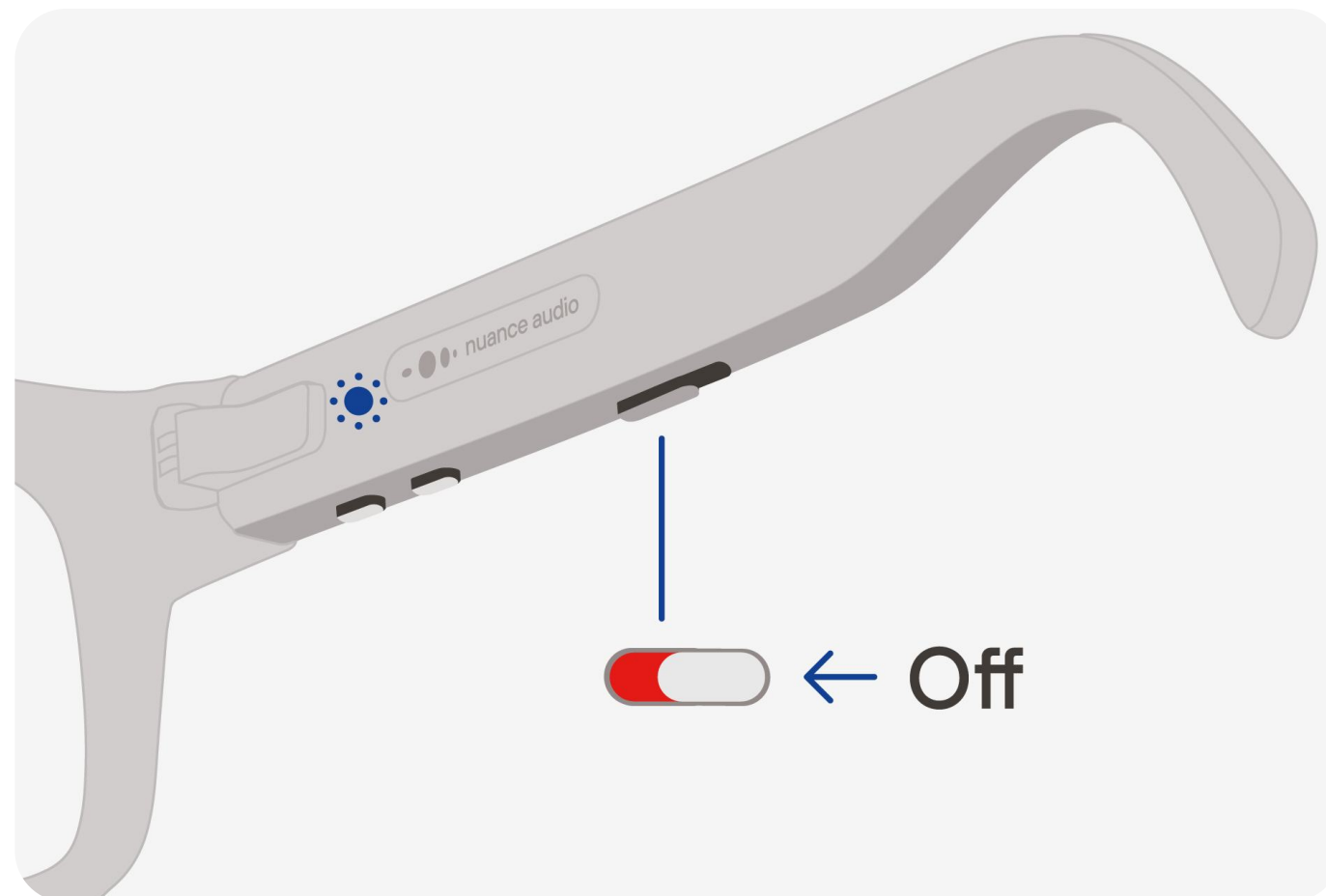

Factory reset

Press (+) and (-) and hold both 10 seconds.

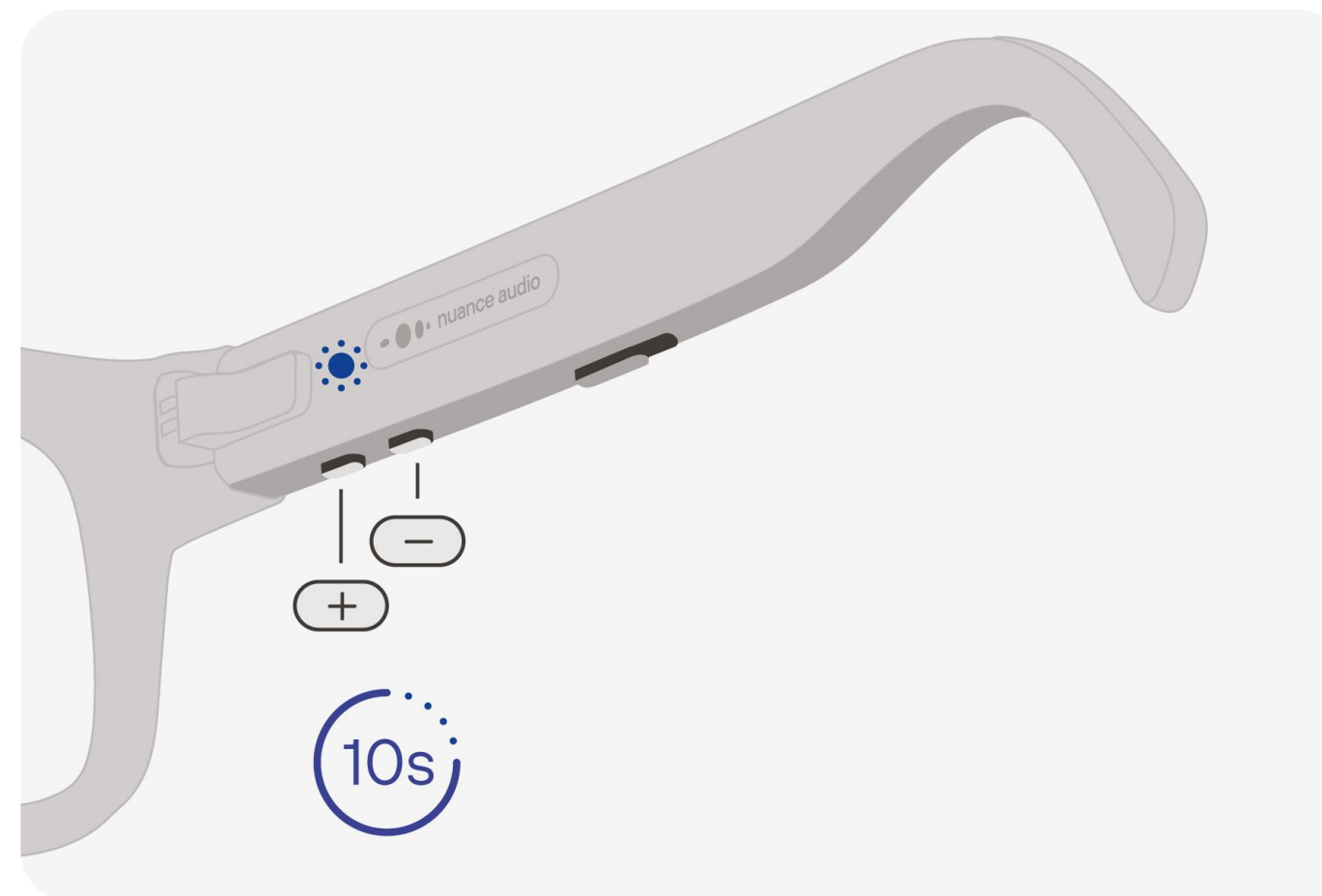



On frame interactions

General

Sleep
Fold



Wake up
Unfold






On frame interactions

Volume

Volume up
Press (+)



Volume down
Press (-)

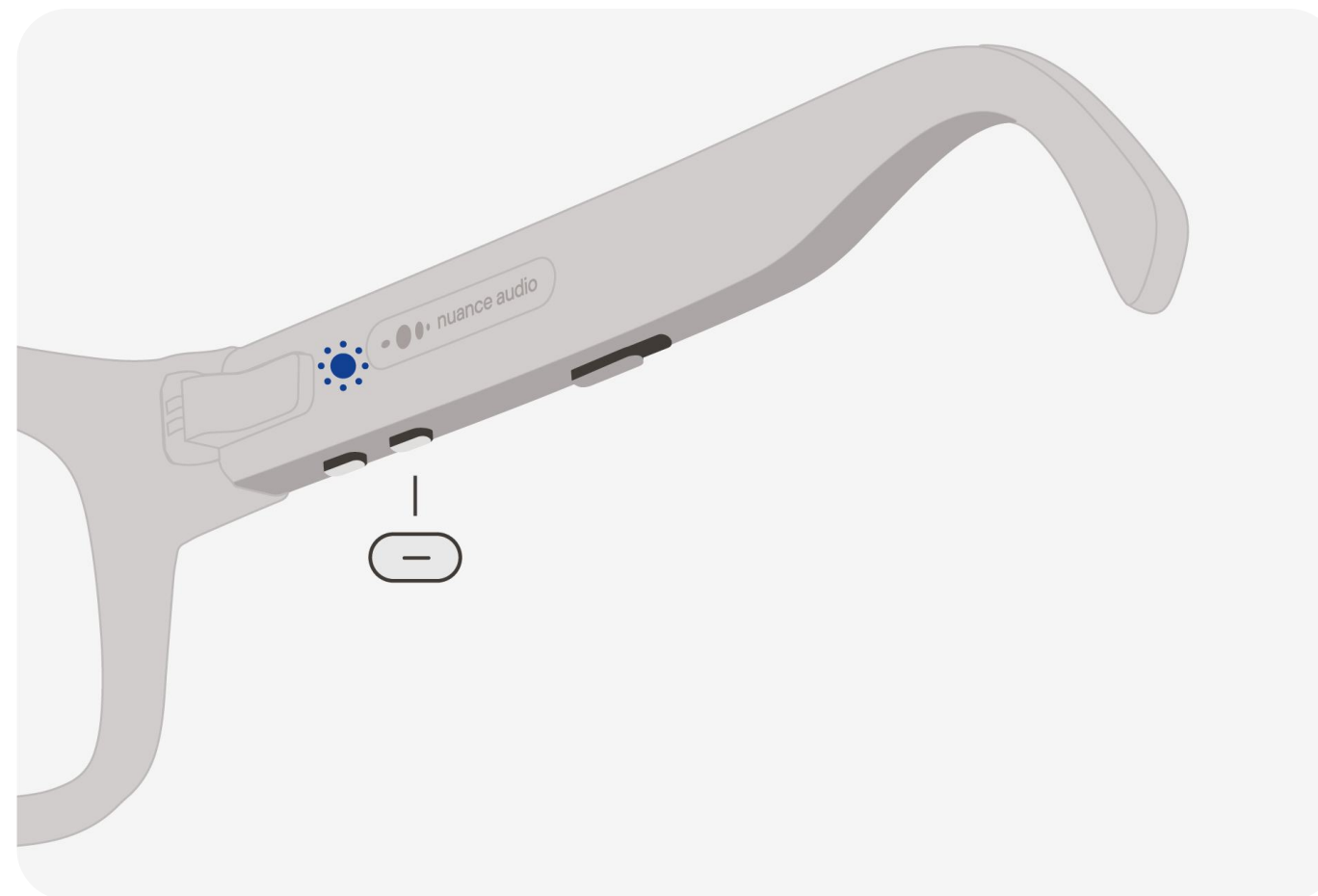



nuance audio

in tune with life

For the most updated Instructions
for Use please follow this link

[↗ View instruction for use](#)